

## Captia™ ENA Combined Screen

REF	2338570	96 Tests
REF	2338571	480 Tests

Pour d'autres langues  
Für andere Sprachen  
Para otras lenguas  
Per le altre lingue  
Dla innych języków

Para outras línguas  
Για τις άλλεςλώσεις  
För andra språk  
For andre språk  
Pro jiné jazyky



[www.trinitybiotech.com](http://www.trinitybiotech.com)

### INTENDED USE

The Trinity Biotech Captia™ ENA Combined Screen Enzyme-Linked Immunosorbent Assay (ELISA) is intended for the detection of antibodies to Sm, RNP, Ro, La, Scl-70, and Jo-1 in human sera. The assay is to be used to detect antibodies in a single serum specimen. The results of the assay are to be used as an aid to the diagnosis of systemic lupus erythematosus, mixed connective tissue disease, Sjogren syndrome, scleroderma, polymyositis, and dermatomyositis. This assay does not differentiate these conditions. **For *in vitro* diagnostic use. High complexity test.**

### INTRODUCTION

Systematic autoimmune disease is characterized by the presence of circulating autoantibodies directed to a wide variety of cellular antigens.<sup>1,2,3</sup> Systemic lupus erythematosus (SLE), commonly referred to as Lupus, is the best known of these diseases. Other possible connective tissue diseases include mixed connective tissue disease (MCTD), Sjogren syndrome, scleroderma, and polymyositis/dermatomyositis. The majority can be diagnosed by clinical presentation and their antibody profiles to the various antigens involved, which include dsDNA, Sm, RNP, Ro, La, Scl-70, Jo-1, and Histones. Therefore, immunoassays for autoantibodies are useful for diagnostic and prognostic evaluations of autoimmune disease.<sup>1,2,3</sup>

The Sm (Smith) antigen is composed of nuclear RNA and several polypeptides. Antibodies to Sm are present in approximately 30% of patients with SLE. Sm is a very specific marker for SLE. Sm antibodies are very rare in other autoimmune diseases and normals.<sup>1,2,3</sup>

RNP antigen is composed of RNA and several protein components. The antigen is closely associated with Sm antigen. Approximately 45% of patients with SLE have antibodies to RNP. RNP antibodies are also found in patients with discoid lupus, Sjogren syndrome, scleroderma, rheumatoid arthritis and MCTD. Very high titers to RNP is a characteristic of MCTD.<sup>1,2,3</sup>

Ro (SS-A) antigen consists of protein complexed with RNA. Approximately 60% of patients with Sjogren syndrome and 35% of SLE patients have antibodies to Ro. SLE patients with Ro antibodies alone (negative La/SSB) have a greater chance of developing nephritis.<sup>1,2,3</sup>

La (SS-B) antigen consists of a protein complexed to RNAs. Approximately 60% of patients with Sjogren syndrome and 15% of SLE patients have antibody to La. Ro (SSA) antibody occurs in virtually all sera with La antibody. However, some SLE patients have Ro only antibody profiles.<sup>1,2,3</sup>

Scl-70 antigen is a nuclear enzyme topoisomerase 1. Approximately 20-30% of patients with scleroderma have antibodies to Scl-70. The antibody is rarely seen in other conditions. Scl-70 antibodies are usually present in the diffuse subtype of scleroderma while antibodies to centromere proteins are dominant in the limited subtype of scleroderma or CREST syndrome.<sup>1,2,3</sup>

Jo-1 antigen is the cellular enzyme histidyl-tRNA synthetase. Approximately 30% of patients with polymyositis (PM) and 10% of patients with dermatomyositis (DM) have antibodies to Jo-1. Antibodies to Jo-1 are seen rarely in normals and other diseases except for those which overlap with PM-DM.<sup>1,2,3</sup>

Classically, antibodies to autoimmune antigens are detected by double immunodiffusion. However, the test is lengthy and suffers weak sensitivity. Enzyme-Linked Immunosorbent Assays (ELISAs) combine greater sensitivity with ease of use. Many ELISAs have been developed and validated for detecting autoantibodies to various antigens.<sup>4,5,6,7</sup>

### PRINCIPLE OF THE ASSAY

The Trinity Biotech ENA Combined Screen test is an Enzyme-Linked Immunosorbent Assay to detect IgG, IgA, and IgM antibodies to ENA antigens. Purified ENA antigens (Sm, RNP, Ro, La, Scl-70, Jo-1) are attached to a solid phase microassay well. Enzyme-Linked Immunosorbent Assays (ELISA) rely on the ability of biological materials (i.e., antigens) to adsorb to plastic surfaces such as polystyrene (solid phase). When antigens bound to the solid phase are brought into contact with a patient's serum, antigen specific antibody, if present, will bind to the antigen on the solid phase forming antigen-antibody complexes. Excess antibody is removed by washing. This is followed by the addition of goat anti-human IgG,A,M conjugated with horseradish peroxidase which then binds to the antibody-antigen complexes. The excess conjugate is removed by washing, followed by the addition of Chromogen/Substrate, tetramethylbenzidine (TMB). If specific antibody to the antigen is present in the patient's serum, a blue color develops. When the enzymatic reaction is stopped with 1N H<sub>2</sub>SO<sub>4</sub>, the contents of the wells turn yellow. The color, which is indicative of the concentration of antibody in the serum, can be read on a suitable spectrophotometer or ELISA microwell plate reader.<sup>10,11,12,13</sup>

### KIT PRESENTATION

#### MATERIALS SUPPLIED

Each kit contains the following components in sufficient quantities to perform the number of tests indicated on the package label.

- Sm, RNP, Ro, La, Scl-70, Jo-1 antigen coated microassay plate : 96 wells, configured in twelve 1x8 strips, stored in a foil pouch with desiccant. (96T: one plate; 480T: five plates)
- Serum Diluent Type III: Ready to use. Contains buffer, BSA and Tween-20, and ProClin® (0.1%) as preservative. (96T: one bottle, 30 mL; 480T: two bottles, 60 mL each)
- High Positive Control: Human serum or defibrinated plasma. Sodium azide (< 0.1%) and pen/strep (0.01%) added as preservatives, with established range printed on vial label. The High Positive Control is utilized to control the upper dynamic range of the assay. (96T: one vial, 0.4 mL; 480T: one vial, 0.8 mL) \*
- Calibrator: Human serum or defibrinated plasma. Sodium azide (< 0.1%) and pen/strep (0.01%) added as preservatives, with kit specific factor printed on vial label. The Calibrator is used to calibrate the assay to account for day-to-day fluctuations in temperature and other testing conditions. (96T: one vial, 0.4 mL; 480T: one vial, 0.8 mL) \*
- Negative Control: Human serum or defibrinated plasma. Sodium azide (< 0.1%) and pen/strep (0.01%) added as preservatives, with established range printed on vial label. The Negative Control is utilized to control the negative range of the assay. (96T: one vial, 0.4 mL; 480T: one vial, 0.8 mL) \*
- Low Positive Control: Human serum or defibrinated plasma. Sodium azide (< 0.1%) and pen/strep (0.01%) added as preservatives, with range printed on vial label. The Low Positive Control is utilized to control the range near the cutoff of the assay. (96T: one vial, 0.4 mL; 480T: one vial, 0.8 mL) \*
- Horseradish-peroxidase (HRP) Conjugate: Ready to use. Goat anti-human IgG, IgA, and IgM containing ProClin® (0.1%) and gentamicin as preservatives. (96T: one bottle, 15 mL; 480T: five bottles, 15 mL each)
- Wash Buffer Type II (20X concentrate): Dilute 1 part concentrate + 19 parts deionized or distilled water. Contains TBS, Tween-80 and ProClin® (0.1%) as a preservative. (96T: one bottle, 50 mL; 480T: one bottle, 250 mL)
- Chromogen/Substrate Solution Type I: Tetramethylbenzidine (TMB), ready to use. The reagent should remain closed when not in use. If allowed to evaporate, a precipitate may form in the reagent wells. (96T: one bottle, 15 mL; 480T: five bottles, 15 mL each)
- Stop Solution: Ready to use, contains a 1N H<sub>2</sub>SO<sub>4</sub> solution. (96T: one bottle, 15 mL; 480T: five bottles, 15 mL each)

\* Note: serum vials may contain excess volume.

The following components are not kit lot # dependent and may be used interchangeably within the Trinity Biotech ELISA Autoimmune Kits: Serum Diluent Type III, Chromogen/Substrate Solution Type I, Wash Buffer Type II, and Stop Solution. Please check that the appropriate Trinity Biotech Reagent Type (Type I, Type II, etc.) is used for the assay.

### ADDITIONAL REQUIREMENTS

- Wash bottle, automated or semi-automated microwell plate washing system.
- Micropipettes, including multichannel, capable of accurately delivering 10-200 µL volumes (less than 3% CV).
- One liter graduated cylinder.
- Paper towels.
- Test tube for serum dilution.
- Reagent reservoirs for multichannel pipettes.
- Pipette tips.
- Distilled or deionized water (dH<sub>2</sub>O), CAP (College of American Pathology) Type 1 or equivalent.<sup>15,16</sup>
- Timer capable of measuring to an accuracy of +/- 1 second (0 - 60 minutes).
- Disposal basins and 0.5% sodium hypochlorite (50 mL bleach in 950 mL dH<sub>2</sub>O).
- Single or dual wavelength microplate reader with 450 nm filter. If dual wavelength is used, set the reference filter to 600-650 nm. Read the Operator's Manual or contact the instrument manufacturer to establish linearity performance specifications of the reader.

Note: Use only clean, dry glassware.

### STORAGE AND STABILITY

- Store unopened kit between 2° and 8° C. The test kit may be used throughout the expiration date of the kit. Refer to the package label for the expiration date.
- Unopened microassay plates must be stored between 2° and 8° C. Unused strips must be immediately resealed in a sealable bag with desiccant, and returned to storage at 2° and 8° C. If the bag is resealed with tape, the wells are stable for 30 days. If the bag is resealed with a heat sealer, the wells are stable until their labeled expiration date.
- Store HRP Conjugate Solution between 2° and 8° C.
- Store the Calibrator, High Positive, Low Positive and Negative Controls between 2° and 8° C
- Store Serum Diluent Type III and 20X Wash Buffer Type II between 2° and 8° C.
- Store the Chromogen/Substrate Solution I between 2° and 8° C. The reagent should remain closed when not in use. If allowed to evaporate, a precipitate may form in the reagent wells.
- Store 1X (diluted) Wash Buffer Type II at room temperature (21° to 25° C) for up to 5 days, or up to 1 week between 2° and 8° C.

Note: If constant storage temperature is maintained, reagents and substrate will be stable for the dating period of the kit. Refer to package label for expiration date. Precautions were taken in the manufacture of this product to protect the reagents from contamination and bacteriostatic agents have been added to the liquid reagents. Care should be exercised to protect the reagents in this kit from contamination. Do not use if evidence of microbial contamination or precipitation is present.

### PRECAUTIONS

- For *in vitro* diagnostic use.
- The human serum components used in the preparation of the Controls and Calibrator in this kit have been tested by an FDA approved method for the presence of antibodies to human immunodeficiency virus 1 & 2 (HIV 1&2), hepatitis C (HCV) as well as hepatitis B surface antigen and found negative. Because no test method can offer complete assurance that HIV, HCV, hepatitis B virus, or other infectious agents are absent, specimens and human-based reagents should be handled as if capable of transmitting infectious agents.
- The Centers for Disease Control & Prevention and the National Institutes of Health recommend that potentially infectious agents be handled at the Biosafety Level 2.<sup>8</sup>
- The components in this kit have been quality control tested as a Master Lot unit. Do not mix components from different lot numbers except Chromogen/Substrate Solution Type I, Stop

Solution, Wash Buffer Type II, and Serum Diluent Type III. Do not mix with components from other manufacturers.

- Do not use reagents beyond the stated expiration date marked on the package label.
- All reagents must be at room temperature (21° to 25°C) before running assay. Remove only the volume of reagents that is needed. **Do not pour reagents back into vials as reagent contamination may occur.**
- Before opening Control and Calibrator vials, tap firmly on the benchtop to ensure that all liquid is at the bottom of the vial.
- Use only distilled or deionized water and clean glassware.
- Do not let wells dry during assay; add reagents immediately after completing wash steps.
- Avoid cross-contamination of reagents. Avoid splashing or generation of aerosols. Wash hands before and after handling reagents. **Cross-contamination of reagents and/or samples could cause erroneous results.**
- If washing steps are performed manually, wells are to be washed three times. Up to five wash cycles may be necessary if a washing manifold or automated equipment is used.
- Sodium azide inhibits Conjugate activity. Clean pipette tips must be used for the Conjugate addition so that sodium azide is not carried over from other reagents.**
- Certain reagents in this kit contain sodium azide for use as a preservative. It has been reported that sodium azide may react with lead and copper in plumbing to form explosive compounds. When disposing, flush drains with water to minimize build-up of metal azide compounds.
- Never pipette by mouth or allow reagents or patient sample to come into contact with skin. Reagents containing ProClin®, sodium azide, and TMB may be irritating. Avoid contact with skin and eyes. In case of contact, immediately flush area with copious amounts of water.
- If a sodium hypochlorite (bleach) solution is being used as a disinfectant, do not expose to work area during actual test procedure because of potential interference with enzyme activity.
- Avoid contact of Stop Solution (1N sulfuric acid) with skin or eyes. If contact occurs, immediately flush area with copious amounts of water.
- Caution:** Liquid waste at acid pH must be neutralized prior to adding sodium hypochlorite (bleach) solution to avoid formation of poisonous gas. Recommend disposing of reacted, stopped plates in biohazard bags. See Precaution 3.
- Do not use Chromogen/Substrate Solution if it has begun to turn blue.
- The concentrations of anti-extractable nuclear antigens in a given specimen determined with assays from different manufacturers can vary due to differences in assay methods and reagent specificity.

The safety data sheet is available upon request.



#### WARNING

Serum Diluent, Conjugate, and Wash Buffer contain 0.1% ProClin 300®, a biocidal preservative that may cause sensitization by skin contact; prolonged or repeated exposure may cause allergic reaction in certain sensitive individuals.

**H317: May cause an allergic skin reaction.**

**P280:** Wear protective gloves / protective clothing / eye protection / face protection.

**P302 + P352:** IF ON SKIN: Wash with plenty of soap and water.

**P333 + P313:** If skin irritation or rash occurs: Get medical advice/ attention.

**P501:** Dispose of contents and container in accordance to local, regional, national and international regulations.

#### WARNING

Serum Diluent and Controls contain < 0.1% sodium azide.

**H302:** Harmful if swallowed

**P264:** Wash thoroughly with plenty of soap and water after handling

**P270:** Do not eat, drink or smoke when using this product

**P301+P312:** IF SWALLOWED: Call a POISON CENTER or doctor/physician if you feel unwell

**P330:** If swallowed, rinse mouth

**P501:** Dispose of contents/container to in accordance to local, regional, national and international regulations.

#### SPECIMEN COLLECTION AND STORAGE

- Handle all blood and serum as if capable of transmitting infectious agents.<sup>8</sup>
- Optimal performance of the Trinity Biotech ELISA kit depends upon the use of fresh serum samples (clear, non-hemolyzed, nonlipemic, non-icteric). A minimum volume of 50 µL is recommended, in case repeat testing is required. Specimens should be collected aseptically by venipuncture.<sup>9</sup> Early separation from the clot prevents hemolysis of serum.
- Store serum between 2° and 8° C if testing will take place within five days. If specimens are to be kept for longer periods, store at -20° C to -70° C in a non-defrosting freezer. Do not use a frostfree freezer because it may allow the specimens to go through freeze-thaw cycles and degrade antibody. Samples that are improperly stored or are subjected to multiple freeze-thaw cycles may yield erroneous results.
- Serum containing visible particulate matter can be spun down utilizing slow speed centrifugation.
- Do not use heat inactivate sera.
- The NCCLS provides recommendations for storing blood specimens (Approved Standard - Procedures for the Handling and Processing of Blood Specimens, H18-A. 1990).<sup>9</sup>

#### METHODS FOR USE

##### PREPARATION FOR THE ASSAY

- All reagents must be removed from refrigeration and allowed to come to room temperature before use (21° to 25° C). Return all reagents to refrigerator promptly after use.
- All samples and controls should be vortexed before use.
- Dilute 50 mL of the 20X Wash Buffer Type II to 1L with distilled and/or deionized H<sub>2</sub>O. Mix well.

##### ASSAY PROCEDURE

- Place the desired number of strips into a microwell frame. Allow six (6) Control/Calibrator determinations (one Negative Control, three Calibrators, one High Positive Control and one Low Positive Control) per run. A reagent blank (RB) should be run on each assay. Check

software and reader requirements for the correct Control/Calibrator configuration. Return unused strips to the sealable bag with desiccant, seal and immediately refrigerate.

#### Example Configuration:

Plate Location	Sample Description	Plate Location	Sample Description
1A	RB	2A	Patient #2
1B	NC	2B	Patient #3
1C	Cal	2C	Patient #4
1D	Cal	2D	Patient #5
1E	Cal	2E	Patient #6
1F	HPC	2F	Patient #7
1G	LPC	2G	Patient #8
1H	Patient #1	2H	Patient #9

RB = Reagent Blank - Well without serum addition run with all reagents. Utilized to blank reader.

NC = Negative Control

Cal = Calibrator

HPC = High Positive Control

LPC = Low Positive Control

- Dilute test sera, Calibrator and Control sera 1:21 (e.g., 10 µL + 200 µL) in Serum Diluent. Mix well. (For manual dilutions it is suggested to dispense the Serum Diluent into the test tube first and then add the patient serum.)
- To individual wells, add 100 µL of the appropriate diluted Calibrator, Controls and patient sera. Add 100 µL of Serum Diluent to reagent blank well. Check software and reader requirements for the correct reagent blank well configuration.
- Incubate each well at room temperature (21° to 25° C) for **30 minutes +/- 1 minute.**
- Aspirate or shake out liquid from all wells. If using semi-automated or automated washing equipment add 250-300 µL of diluted Wash Buffer to each well. Aspirate or shake to remove all liquid. Repeat the wash procedure two times (for a total of three (3) washes) for manual or semi-automated equipment or four times (for a total of five (5) washes) for automated equipment. After the final wash, blot the plate on paper toweling to remove all liquid from the wells.

**\*\*IMPORTANT NOTE:** Regarding steps 5 and 8 - Insufficient or excessive washing will result in assay variation and will affect validity of results. Therefore, for best results the use of semi-automated or automated equipment set to deliver a volume to completely fill each well (250-300 µL) is recommended. A total of up to five (5) washes may be necessary with automated equipment. **Complete removal of the Wash Buffer after the last wash is critical for the accurate performance of the test. Also, visually ensure that no bubbles are remaining in the wells.**

- Add 100 µL Conjugate to each well, including reagent blank well. Avoid bubbles upon addition as they may yield erroneous results.
- Incubate each well at room temperature (21° to 25° C) for **30 minutes +/- 1 minute.**
- Repeat wash as described in Step 5.
- Add 100 µL Chromogen/Substrate Solution (TMB) to each well, including the reagent blank well, maintaining a constant rate of addition across the plate.
- Incubate each well at room temperature (21° to 25° C) for **15 minutes +/- 1 minute.**
- Stop reaction by addition of 100 µL of Stop Solution (1N H<sub>2</sub>SO<sub>4</sub>) following the same order of Chromogen/Substrate addition, including the reagent blank well. Tap the plate gently along the outsides, to mix contents of the wells. The plate may be held up to 1 hour after addition of the Stop Solution before reading.
- The developed color should be read on an ELISA plate reader equipped with a 450 nm filter. If dual wavelength is used, set the reference filter to 600-650 nm. The instrument should be blanked on air. The reagent blank must be less than 0.150 Absorbance at 450 nm. If the reagent blank is > 0.150 the run must be repeated. Blank the reader on the reagent blank well and then continue to read the entire plate. Dispose of used plates after readings have been obtained.

#### QUALITY CONTROL

For the assay to be considered valid the following conditions must be met:

- Calibrators and Controls must be run with each test run.
- Reagent Blank must be < 0.150 O.D. (Optical Density) at 450 nm (when read against Air Blank).
- The mean O.D. value for the Calibrator should be ≥ 0.300 at 450 nm (when read against Reagent Blank).
- The Index Values for the High Positive, Low Positive, and Negative Controls should be in their respective ranges printed on the vials. If the control values are not within their respective ranges, the test should be considered invalid and should be repeated.
- Additional Controls may be tested according to guidelines or requirements of local, state, and/or federal regulations or accrediting organizations.
- Refer to NCCLS C24A for guidance on appropriate Quality Control practices.<sup>14</sup>
- If above criteria are not met on repeat, contact Trinity Biotech Technical Service.

#### INTERPRETATION

- Mean Calibrator O.D. - Calculate the mean value for the Calibrator from the three Calibrator determinations. If any of the three Calibrator Values differ by more than 15% from the mean, discard that value and calculate the average of the two remaining values.
- Correction Factor - To account for day-to-day fluctuations in assay activity due to room temperature and timing, a Correction Factor is determined by Trinity Biotech for each lot of kits. The Correction Factor is printed on the Calibrator vial.
- Cutoff Calibrator Value - The Cutoff Calibrator Value for each assay is determined by multiplying the Correction Factor by the mean Calibrator O.D. determined in Step 1.
- Index Value - Calculate an Index Value for each specimen by dividing the specimen O.D. value by the Cutoff Calibrator Value determined in Step 3.

**Example:**

O.D.s obtained for Calibrator	= 0.38, 0.42, 0.40
Mean O.D. for Calibrator	= 0.40
Correction Factor	= 0.50
Cutoff Calibrator Value	= 0.50 x 0.40 = 0.20
O.D. obtained for patient sera	= 0.60
Index Value	= 0.60/0.20 = 3.00

**ANALYSIS**

1. The patients' Index Values are interpreted as follows:

Index Value	Results	Interpretation
≤ 0.90	Negative	No detectable antibody to Extractable Nuclear Antigens by the ELISA test.
0.91-1.09	Equivocal	Samples should be retested. See number 2 below.
≥ 1.10	Positive	Indicates presence of detectable antibody to Extractable Nuclear Antigens by the ELISA test.

2. Samples that remain equivocal after repeat testing should be retested on an alternate method or test a new sample.

**EXPECTED VALUES**

- From 1-5% of apparently normal individuals may contain autoantibodies.<sup>12,3</sup>
- Antibodies to the various ENA's vary depending on the disease state. The following are the reported incidences of ENA antibodies in various diseases.

Antibody	DISEASE					
	SLE	SS	PSS	MCTD	PM	DM
Sm	~ 30%	-	-	~ 2%	-	-
RNP	~ 45%	~ 5%	~ 25%	~ 95%	-	-
Ro	~ 30%	~ 60%	-	-	-	-
La	~ 15%	~ 60%	-	-	-	-
Scl-70	-	-	~ 25%	-	-	-
Jo-1	-	-	-	-	~ 30%	~ 10%

SLE - systemic lupus erythematosus  
 SS - Sjogren syndrome  
 PSS - scleroderma  
 MCTD - mixed connective tissue disease  
 PM - polymyositis  
 DM - dermatomyositis

**LIMITATIONS OF USE**

- The result of the assay should not be interpreted as being diagnostic. The results should only be used as an aid to diagnosis. The results should be interpreted in conjunction with the clinical evaluation of the patient.
- Sera from patients with other autoimmune diseases and from normal individuals may contain autoantibodies.
- The assay should be used only with serum. Icteric, lipemic, hemolyzed and heat inactivated serum should be avoided.
- Index Values of ≥ 10.00 should be reported as greater than 10.

**PERFORMANCE CHARACTERISTICS**

**SENSITIVITY AND SPECIFICITY**

The Trinity Biotech ENA Combined Screen ELISA kit was evaluated relative to a commercially available ELISA test kit. Table 1 summarizes the data.

**Table 1: Sensitivity and Specificity of the Trinity Biotech ENA Combined Screen ELISA Kit**

		Trinity Biotech ENA Combined Screen ELISA Kit			Total
		Positive ≥ 1.10	Equival 0.91-1.09	Negative ≤ 0.90	
Alternate ELISA Kit	Positive > 25	70	-	-	70
	Equivocal 21-25	1	-	-	1
	Negative < 21	3	-	128	131
	Total	74	0	128	202

Sera falling in the equivocal range were excluded from the following calculations.

Relative Sensitivity = 70/73 = 96 %  
 Relative Specificity = 128/128 = 100 %  
 Relative Agreement = 198/201 = 98.5 %

\* The positive sera consisted of the following specificities Sm-18, RNP- 17, Ro-11, La-6, Scl-70-10, Jo-1-12

**REPRODUCIBILITY**

The precision of the Trinity Biotech ENA Combined Screen ELISA kit was determined by testing seven different sera eight times each on three different assays. The data are summarized in Table 2. With proper technique, the user should obtain C.V.'s of less than 15%.

**Table 2: Precision Data ENA Combined Screen ELISA Assay**

Serum#	Assay 1 (n=8)			Assay 2 (n=8)			Assay 3 (n=8)			Inter Assay (n=24)		
	X	S.D.	C.V.	X	S.D.	C.V.	X	S.D.	C.V.	X	S.D.	C.V.
1	6.89	0.319	4.63%	7.29	0.351	4.81%	7.46	0.247	3.31%	7.21	0.38	5.31%
2	3.24	0.416	12.8%	3.40	0.251	7.38%	3.67	0.196	5.34%	3.44	0.342	9.94%
3	6.24	0.253	4.05%	6.83	0.249	3.65%	6.94	0.307	4.42%	6.67	0.407	6.10%
4	5.67	0.186	3.28%	5.89	0.161	2.74%	6.24	0.271	4.34%	5.93	0.312	5.26%
5	5.95	0.382	6.42%	6.04	0.321	5.31%	6.36	0.202	3.18%	6.12	0.348	5.69%
6	0.17	0.065	38.2%	.273	0.061	22.3%	.197	0.098	33.7%	.261	0.071	27.2%
7	0.22	0.031	14.2%	.260	0.195	75.0%	.291	0.170	86.3%	.206	0.149	72.3%

X = Mean ENA Value  
 S.D. = Standard Deviation  
 C.V. = Coefficient of Variation

**LINEARITY**

The Trinity Biotech ENA Combined Screen Index Values were determined for serial twofold dilutions of five positive sera. The Index Values were compared to log<sup>2</sup> of dilution by standard linear regression. The data in Table 3 indicates that the assay is linear.

**Table 3: Linearity for the Trinity Biotech ENA Combined Screen ELISA Test**


Serum #	Neat	1:2	1:4	1:8	1:16	1:32	1:64	1:128	r <sup>2</sup>
1	6.53	5.59	4.61	3.78	3.10	2.25	1.38	0.78	0.997
2	6.40	4.83	3.45	2.38	1.57	1.07	0.74		0.944
3	5.21	4.23	3.53	2.94	2.40	1.93	1.43	1.01	0.982
4	8.97	8.22	7.39	5.98	4.61	3.14	1.87	1.14	0.991
5	8.05	7.44	6.53	5.24	4.17	2.52	1.68	1.09	0.988

r<sup>2</sup> = coefficient of determination.

**REFERENCES**

- Tan, E.M. 1982. Autoantibodies to Nuclear Antigens (ANA): Their Immunobiology and Medicine. *Adv. Immunol.* 33: 167-240.
- Nakamura, R.M., and E.M. Tan. 1986. Recent Advances in Laboratory Tests and the Significance of Autoantibodies to Nuclear Antigens in Systemic Rheumatic Disease. *Clin. Lab Med.* 6: 41-53.
- McCarty, G.A., D.W. Valencia, and M.J. Fritzler. 1984. Antinuclear Antibodies: Contemporary Techniques and Clinical Application to Connective Tissue Disease. In: *Antinuclear Antibodies*. Oxford Univ. Press, New York. pp 1-95.
- Kurki, P., M. Gripenberg, P. Parlanen, and T. Helve. 1987. Screening Test for Rheumatic Disease: A Combined Enzyme Immunoassay of Rheumatoid Factors and Autoantibodies to DNA and Extractable Nuclear Antigens. *J. Clin. Pathol.* 40 : 1475-1480.
- Geisler, C. and M. Hoier-Madsen. 1985. An Enzyme-Linked Immunosorbent Assay for Autoantibodies Against the Nuclear Protein Scl-70. *J. Immuno. Methods.* 80: 211-219.
- Biswas, T., F.W. Miller, Y. Takagaki, and P.H. Poltz. 1987. An Enzyme-Linked Immunosorbent Assay for the Detection and Quantitation of Anti-Jo-1 Antibody in Human Serum. *J. Immun. Meth.* 98: 243-248.
- Rubin, R.L. 1992. Enzyme-Linked Immunosorbent Assay for Anti- DNA and Antihistone Antibodies Including Anti-(H2A-H2B). In: *Manual of Clinical Laboratory Immunology*. N.R. Rose, et al., eds. 4th Edition. 108: 735-740.
- CDC-NIH Manual. 1993. In: *Biosafety in Microbiological and Biomedical Laboratories*, 3rd Edition. U. S. Dept. of Health and Human Services, Public Health Service. pp 9-12.
- National Committee for Clinical Laboratory Standards. 1990. Procedures for the Collection of Diagnostic Blood Specimens by Venipuncture Approved Standard. NCCLS Publication H18-A.
- Engvall, E., and P. Perlman. 1971. Enzyme-Linked Immunosorbent Assay, (ELISA) Quantitative Assay of Immunoglobulin G. *Immunochemistry.* 8: 871-874.
- Engvall, E., and P. Perlman. 1971. Enzyme-Linked Immunosorbent Assay, ELISA. Peeters, H., ed. In: *Protides of the Biological Fluids*. Proceedings of the Nineteenth Colloquium, Brugge, Oxford. Pergamon Press. pp 553-556.
- Engvall, E., K. Jonsson, and P. Perlman. 1971. Enzyme-Linked Immunosorbent Assay. II. Quantitative Assay of Protein Antigen, Immunoglobulin-G, By Means of Enzyme-Labelled Antigen and Antibody-Coated tubes. *Biochem. Biophys. Acta.* 251: 427-434.
- Van Weeman, B. K. and A.H.W.M. Schuurs. 1971. Immunoassay Using Antigen-Enzyme Conjugates. *FEBS Letter.* 15: 232-235.
- NCCLS. 1991. National Committee for Clinical Laboratory Standard. Internal Quality Control Testing: Principles & Definition. NCCLS Publication C24- A.
- <http://www.cap.org/html/ftpdirectory/checklistftp.html>. 1998. Laboratory General - CAP (College of American Pathology) Checklist (April 1998). pp 28-32.
- NCCLS. 1997. National Committee for Clinical Laboratory Standard. Preparation and Testing of Reagent Water in the Clinical Laboratory. NCCLS Publication C3- A3.

CAPTIA™ ENA COMBINED SUMMARY OF ASSAY PROCEDURE

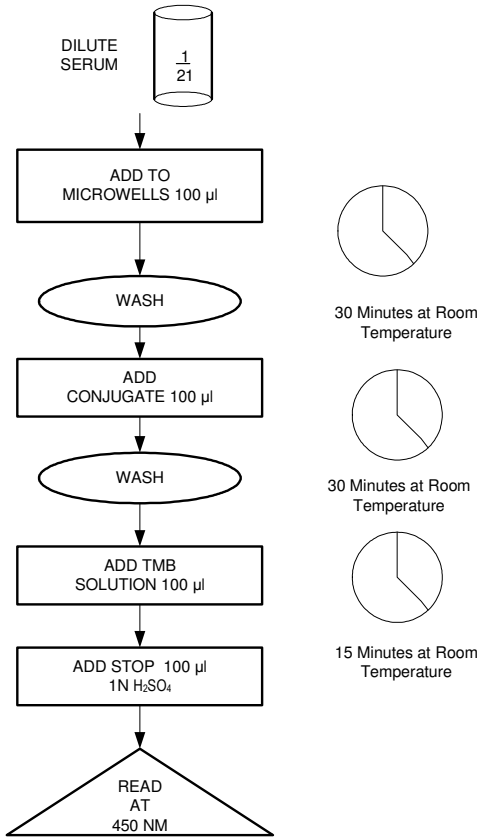

The safety data sheet is available upon request.



**WARNING**

Serum Diluent, Conjugate, and Wash Buffer contain 0.1% ProClin 300®, a biocidal preservative that may cause sensitization by skin contact; prolonged or repeated exposure may cause allergic reaction in certain sensitive individuals.

**H317: May cause an allergic skin reaction.**

**P280: Wear protective gloves / protective clothing / eye protection / face protection.**

**P302 + P352: IF ON SKIN: Wash with plenty of soap and water.**

**P333 + P313: If skin irritation or rash occurs: Get medical advice/ attention.**

**P501: Dispose of contents and container in accordance to local, regional, national and international regulations.**

**WARNING**

Serum Diluent and Controls contain < 0.1% sodium azide.

**H302: Harmful if swallowed**

**P264: Wash thoroughly with plenty of soap and water after handling**

**P270: Do not eat, drink or smoke when using this product**

**P301+P312: IF SWALLOWED: Call a POISON CENTER or doctor/physician if you feel unwell**

**P330: If swallowed, rinse mouth**

**P501: Dispose of contents/container to in accordance to local, regional, national and international regulations.**

ORDERING INFORMATION

KIT		Captia™ ENA Combined Screen Test Kit	
Catalog No.	Item	Quantity	
2338570	Captia™ ENA Combined Screen Test Kit	96 Tests	
2338571	Captia™ ENA Combined Screen Test Kit	480 Tests	

Manufactured	High Pos or Positive Control
Authorized Representative	Low Pos or Cut-Off Control
Consult accompanying documents	Negative Control
Product Number	Calibrator
Lot	Coefficient Factor
Use by	Range
Caution, consult accompanying documents	Standard
Store at 2-8°C	For In Vitro Diagnostic use
Store at 2-30°C	or  Warning
Trinity Biotech USA Jamestown, NY 14701 Tel. 1 800-325-3424 Fax: 716-488-1990	Trinity Biotech plc Bray Co. Wicklow, Ireland Tel. 353 1 2769800 Fax 353 1 2769888 www.trinitybiotech.com

8570-29 Rev L  
03/2015